

OpenFlexure Microscope

Abstract: A 3D-printed microscope, including a flexure-based XYZ stage with sub-micron precision and 8mm range.

Description:

The OpenFlexure microscope is a customisable optical microscope, using either very cheap webcam optics or lab quality, RMS standard microscope objectives. It uses an inverted geometry, and has a high quality mechanical stage which can be motorised using low cost geared stepper motors.

Optomechanics is a crucial part of any microscope; when working at high magnification, it is absolutely crucial to keep the sample steady and to be able to bring it into focus precisely. Accurate motion control is extremely difficult using printed mechanical parts, as good linear motion typically requires tight tolerances and a smooth surface finish. This design for a 3D printed microscope stage uses plastic flexures, meaning its motion is free from friction and vibration. It achieves steps well below 100nm when driven with miniature stepper motors, and is stable to within a few microns over several days.

This design aims to minimise both the amount of post-print assembly required, and the number of non-printed parts required - partly to make it as easy as possible to print, and partly to maximise stability; most of the microscope (including all the parts with flexures) prints as a single piece. The majority of the expense is in the Raspberry Pi and its camera module; the design requires only around 100g of plastic and a few nuts, bolts and other parts. The optics module (containing the camera and lens) can be easily swapped out or modified, for example to change the magnification/resolution by using a microscope objective, or adding a filter cube for fluorescence (build instructions not done yet).

This design is (c) Richard Bowman 2017, released under the CERN open hardware license. You can also buy a kit or an assembled microscope from [WaterScope](#).

A 3D-printed microscope, including a flexure-based XYZ stage with sub-micron precision and 8mm range.

License: CERN Open Hardware License

Notes:

This is the development version of the microscope; ongoing development lives on [Github](#) and there is an open-access paper in [Review of Scientific Instruments](#).

This microscope is available through [WaterScope](#), as a kit (with or without printed parts and Raspberry Pi) or fully assembled.

Requires:

- **Optics module:** [Basic optics module x 1](#), [High-resolution Optics Module x 1](#)
- **Illumination:** [Illumination Arm \(fixed\) x 1](#), [Adjustable Illumination Arm x 1](#), [Illumination Arm with Condenser x 1](#)
- **Fine positioning control:** [Actuator assembly x 3](#)
- **Mechanical stage:** [Microscope Body x 1](#)
- **Sample clip attachment screws:** [8mm M3 screws, Qty: 3 x 1](#)
- **Sample clips:** [Sample Clips x 1](#)
- **Control Computer:** [Raspberry Pi x 1](#)
- **Printer test object:** [Microscope leg test object x 1](#)

Instruction:

Step 1. First, you will need to print or obtain the 3D printed parts. If you bought a kit, skip this step!

The microscope is designed to print without support material. This is quite important; if you use support material it will require a lot of cleaning up, and you may well damage the mechanism. I usually print with a layer size of 0.24mm on my Ormerod, which takes 8 hours for the main structure. "low" quality on an Ultimaker 2 (0.15mm layers) produced similar results in about 5 hours.

If your printer has a standard-sized bed (180mmx180mm should be fine) then it should be possible to print the complete microscope in one go. I do this if I'm using a machine that is well calibrated and reliable. However, I find that it's often more reliable to print in batches (as small parts at the edge of the print bed can detach and cause it to fail). I would recommend:

- Batch 1: Microscope, illumination and Optics module (main part)
- Batch 2: Feet, gears, camera cover, camera board gripper, camera lens remover, gear riser

There is a test file that prints a single leg of the microscope - the Microscope leg test object. It's worth printing this first to check your settings are OK.

Step 2. Make sure you have all the necessary parts and tools. The parts should all be listed in the bill of materials, but for clarity you will need:

Plastic parts:

- 1x microscope body
- 1x optics module
- 1x illumination and rear foot
- 2x tilted foot
- 1x untilted foot
- 3x gear
- 1x camera cover
- 2x sample clip (optional)



Plastic tools (optional but helpful):

- 1x camera board gripper (for v1 or v2 of the camera board)
- 1x camera lens remover (for v1 or v2)
- 1x gear elevator (optional)

Metal hardware:

- 3x M3 hexagon head 25mm screws, stainless steel
- 3x M3 brass nut
- 3x M3 stainless steel washer
- 2x M3 8mm cap head screw (optional, for sample clips)
- 2x M2 6mm cap head screw (optional, to secure camera more strongly)

Electronic parts:

- White LED, resistor, wire, and 2-way JST header connector (assembled as one cable in the kit)
- Raspberry Pi camera module (v2, though v1 works if you substitute the relevant STL files)
- Raspberry Pi (with associated power supply, keyboard, monitor, etc.)

Tools (not supplied in kit):

- 2.5mm hex key (optional, for attaching sample clips)
- 1.5mm hex key (optional, for double-securing camera board)
- tape (electrical tape or PTFE plumbers tape work, though regular sticky tape is also fine)
- sharp craft knife (for trimming tape)
- 3mm drill bit in hand chuck (if you printed the parts yourself and need to open out the holes)

Don't forget the raspberry pi, camera module, and associated screen, power supply, SD card, keyboard, mouse, etc. (I have not listed these explicitly, but they're needed to run the Pi).


Step 3. If you printed the parts yourself, start by opening out the three holes in the microscope body with a drill as shown. Make sure to go all the way through. If you don't have a drill, you can improvise by screwing in an M3 screw all the way, then forcibly rotating it with a screwdriver. While you're at it, do the same for the illumination mount. Also, remove any loose strings of plastic from the underside of the actuator housings or sample stage, using a pair of pliers. The last step shouldn't be necessary if your machine is calibrated nicely for printing bridges. If you purchased a kit, this may well have been done for you.



Step 4. Depending on which version of the microscope body you have, there may be thin supports linking the mount for the optics module to the frame. Snap these off with a craft knife, screwdriver, or elastic band tool.




Step 5. First, assemble the three gears and elastic bands - it's the same procedure done three times, once for each gear. You might find the central column is the best one to start with. NB the instructions for this are in the actuator assembly brick

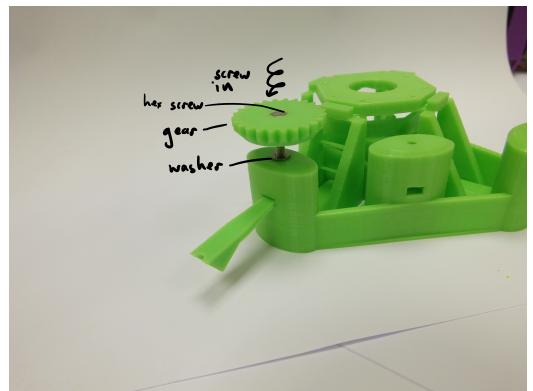


Step 6. Assemble the illumination arm. There are different versions - a fixed one for a single LED and an adjustable one which may have a condenser lens. As with the actuators, the assembly instructions for this are in later sections.





Step 7. Assemble the Raspberry Pi camera into the optics module. There are two versions of this, a high-resolution module using an RMS microscope objective, and the basic optics module described in later sections.

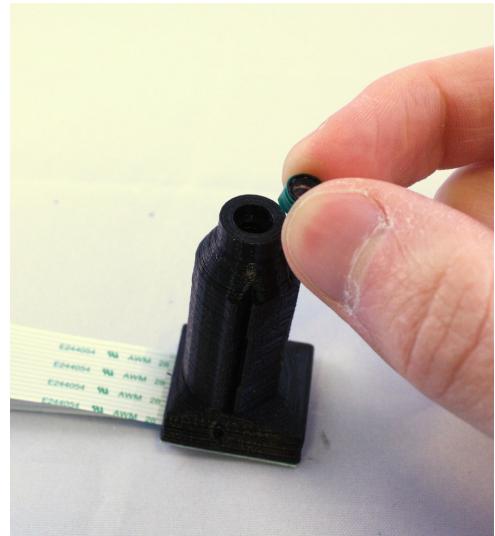

Step 8. Slot the optics module onto the dovetail mount, from the bottom. You may find it easier to fit if you wiggle the module from side to side as it goes in. It may also help to support the dovetail from the top (through the hole in the stage), as it will flex when you push the objective on to it. If you push too hard it is possible to break the flexure mechanism that holds the optics module.



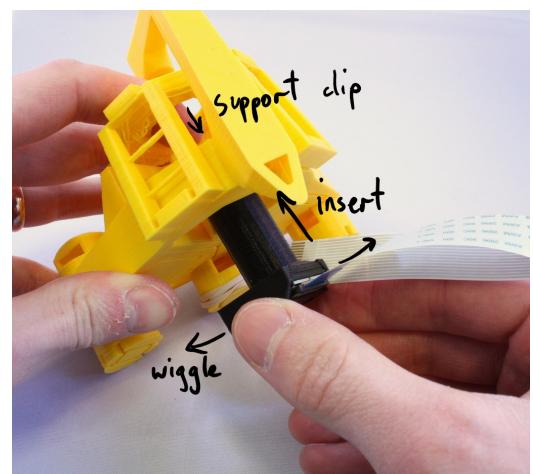


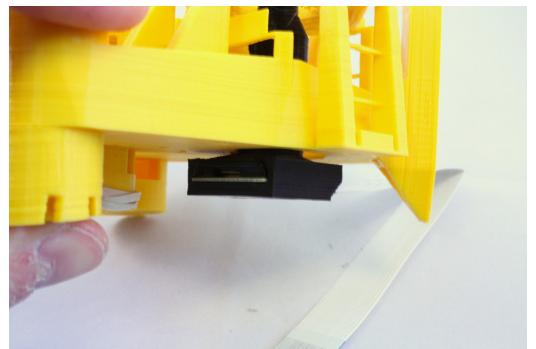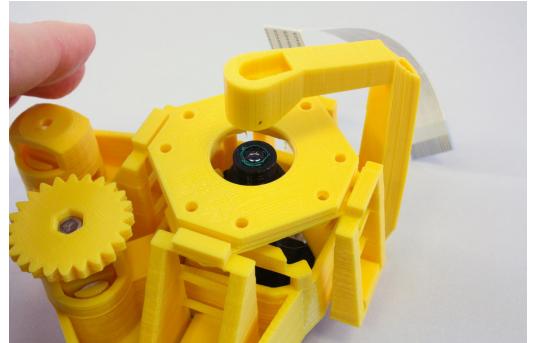


Step 9. Plug the LED cable into the GPIO connector on the Raspberry Pi, to the 0v and 5v lines. These are the second and third pins from the top of the connector, on the outside edge - pins number 4 and 6. Plug in the camera to the camera connector as described in the Raspberry Pi documentation (the connector is next to the Ethernet port, and the contacts on the cable face the port, i.e. they face away from the tab on the plug).





Step 10. After this, there are only the sample clips to go – exactly where you place these will depend on the samples you intend to use, but in any case you simply push the M3 screws into the clips, then screw down into the holes on the stage.

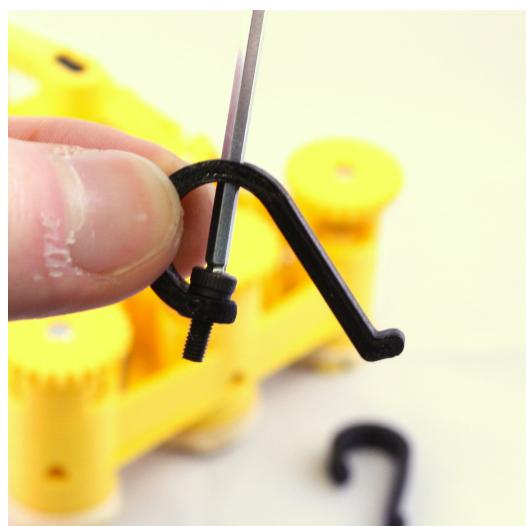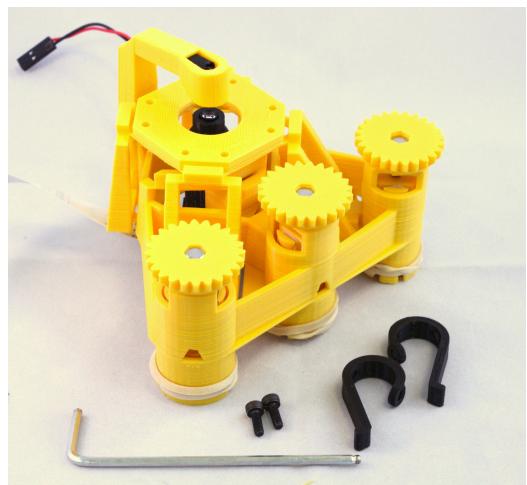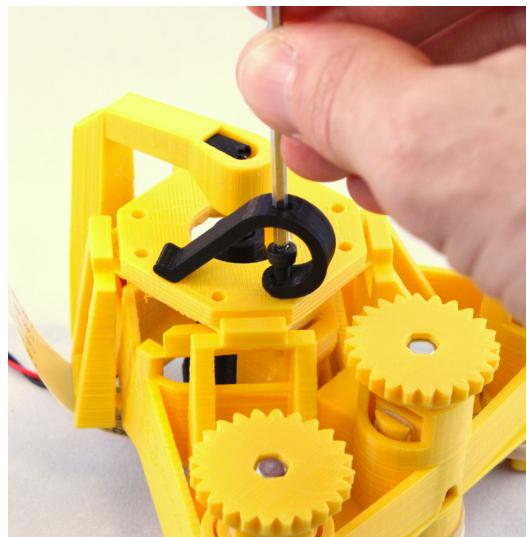

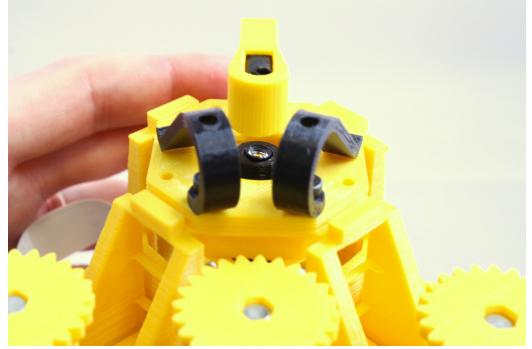

Step 11. Your microscope is now complete – happy observing! You might want to consult the [camera module documentation](#) or [raspicam documentation](#) if you need a hand setting up the camera.

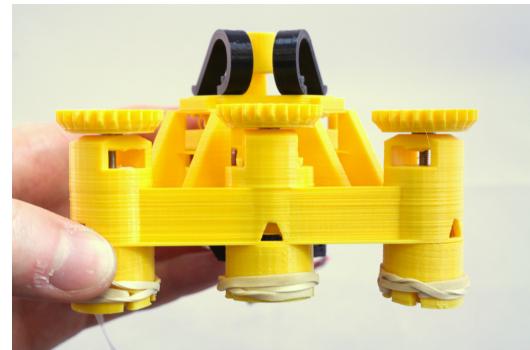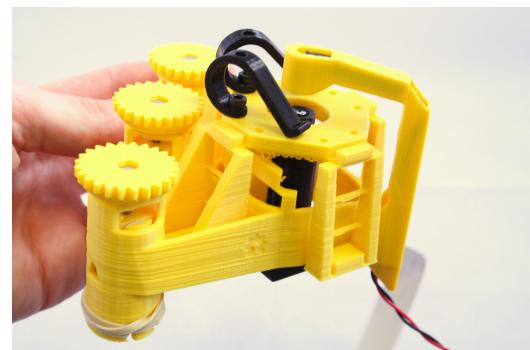

Basic optics module

Abstract: This is a push-fit extension tube that turns a Raspberry Pi camera module into a microscope with a field of view about 400um across and a resolution of around 2um (better than 2um in the case of v2 of the camera module).

Description: This is the optics module for the OpenFlexure Microscope. It holds a Raspberry Pi camera module and lens, but the lens is reversed and placed 40mm from the sensor, allowing it to function as a high-magnification microscope with a resolution of around 2um.

This is a push-fit extension tube that turns a Raspberry Pi camera module into a microscope with a field of view about 400um across and a resolution of around 2um (better than 2um in the case of v2 of the camera module).

Requires:

- **Mechanical holder:** [Optics module plastic parts x 1](#)
- **Optics and Sensor:**
- **Tools for removing the lens:** [Camera lens removal tools x 1](#)

Instruction:

Step 1. First, assemble the necessary tools and parts: the optics module plastic parts (the "objective" and camera board cover), the Raspberry Pi camera module (version 2), the tools for removing the camera lens, optionally two M2 screws to secure the camera. Depending on print quality, you might also need a sharp craft knife or some tape. NB the lens removal tool, board gripper, and optics module are all specific to the camera you're using (v1 or v2 of the Raspberry Pi camera module). This version of the instructions is for version 2 of the camera board, if you have v1 it's a similar process but the older versions of the instructions give more detail. Later production runs of the Raspberry Pi camera module ship with a lens removal tool included, which is better than the printed tool.




Step 2. WARNING! The camera board is static sensitive. Take the usual anti-static precautions (ideally use an anti-static wristband connected to ground, but at the very least make sure you touch an earthed object, such as a metal pipe, before working on the camera module).





Step 3. We need to remove the lens from the camera. To do this, you need the two plastic tools (the board gripper and the lens remover) as well as the camera module. If you are using a USB camera, the lens is held on by two screws from the back, and can be removed by unscrewing these. **SAVE THE SCREWS**, as you will need them later!





Step 4. Remove the protective film from the camera lens.

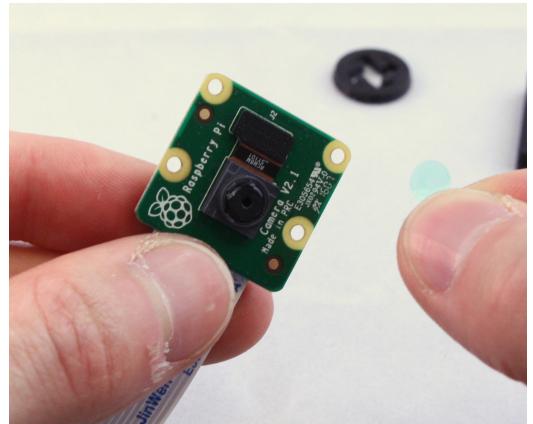

Step 5. There is a small ribbon cable connecting the camera to the PCB that is very easy to break. There is a square plastic jig that fits over the camera and PCB (the "camera board gripper"), which stops the camera twisting and damaging the ribbon cable. Fit this over the camera as shown. Note that the part for v2 of the camera board will sort-of fit v1, but you need to be a little more careful as it's not a perfect fit.

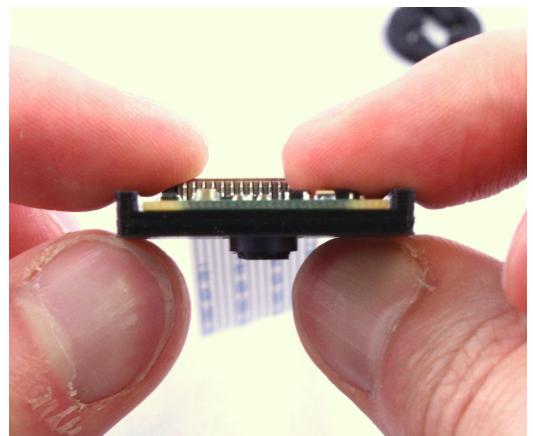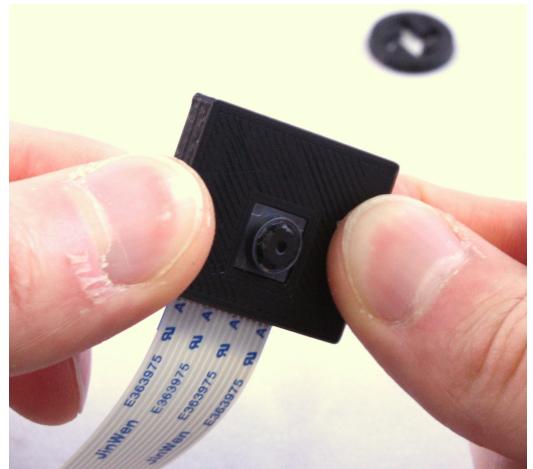




Step 6. Next, unscrew the lens from the camera module. Use the plastic tool to grip the lens module. This is a small circular part with four prongs that fits over the lens of the camera board (version 2 only) as shown. To remove the lens, push the removal tool onto the lens (just the top part, with the little plastic flanges) and turn anticlockwise to remove it. It only works if the prongs are pointing anticlockwise, so make sure it's the right way round. You will feel a crack as it starts to move - this is normal. It's important to use the board gripper to hold the camera chip in place and prevent damage to the delicate ribbon cable. After you've removed the lens, check that the little black ribbon cable connecting the camera module (square plastic housing) to the PCB is still connected - pop it back in by pushing it with a finger if needed. Once you've removed the lens, be sure to place the camera face down on the desk, or put a piece of tape over the square black lens holder; this will help stop dust settling on the sensor, which is hard to clean.

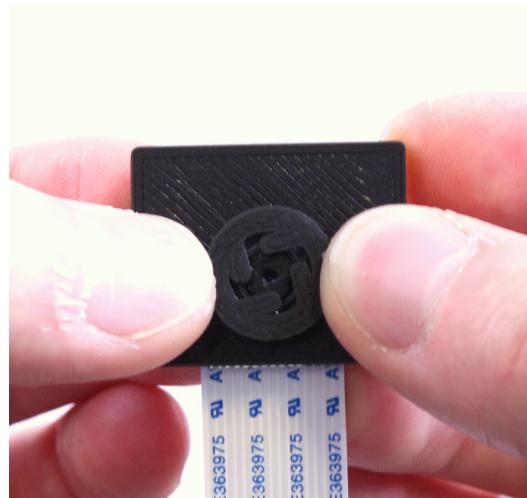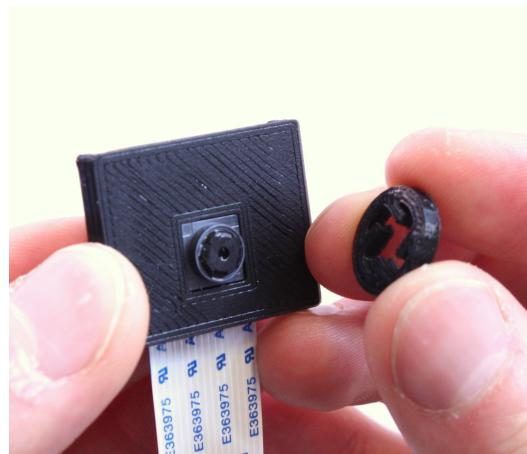



Step 7. Before assembling the parts into the holder, make sure it's free from dust by blowing some air through it, and check there are no strings of plastic in the central hole through the mount.


Step 8. Next, put the lens into the holder. This should just push-fit, but may take a small amount of force, or require a layer of tape wrapped around the lens to make it fit tightly (depending on your printer). If it looks like it will require too much force, you might have a different design of camera module and you may need to re-print the part slightly larger (or get a replacement if you bought a kit). If you wrapped tape around the lens, trim off any tape that protrudes above the lens with a scalpel or sharp craft knife.

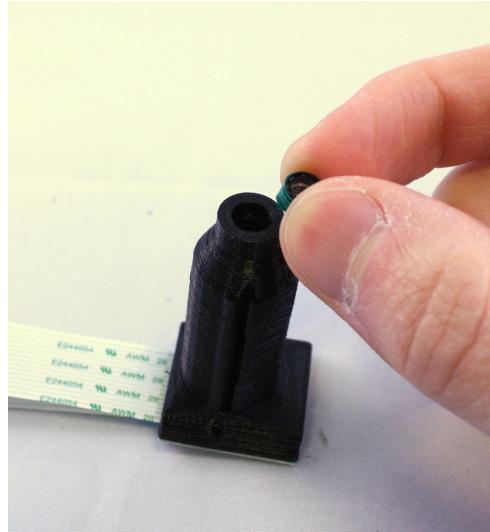

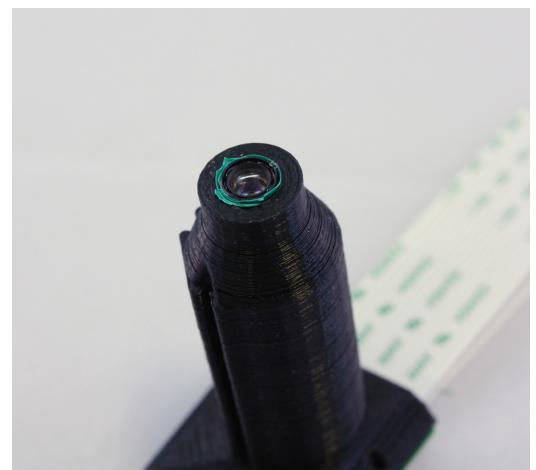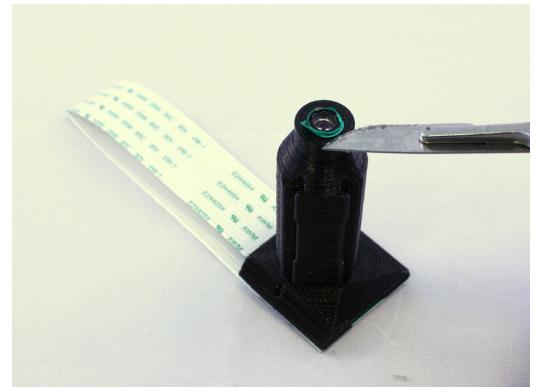


Step 9. Finally, press the camera module onto the bottom of the optics mount. Don't worry if it's not gripped very strongly, the cover will hold it in place.

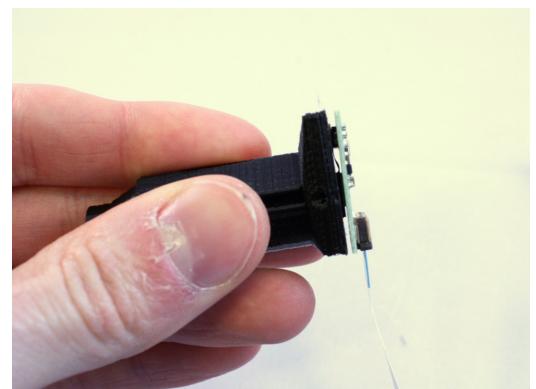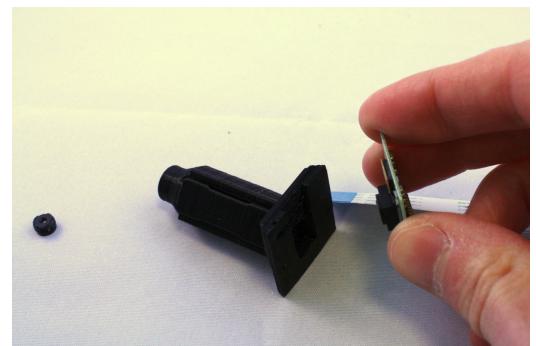

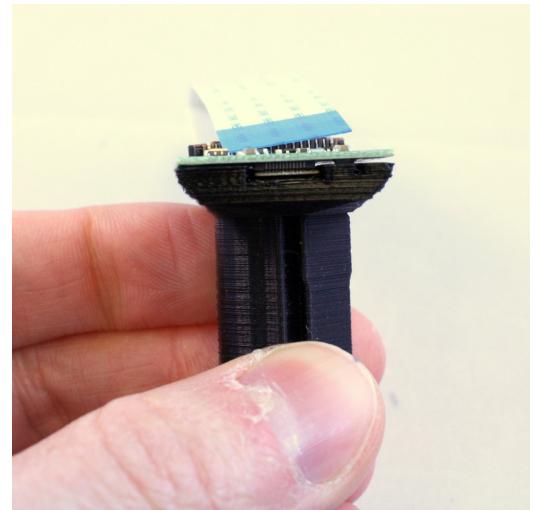

Step 10. Place the cover on top of the camera board, then screw the cover and the camera board onto the optics module, using two M2 screws.

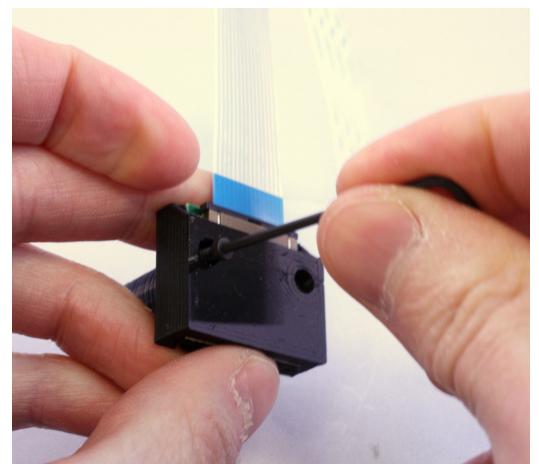

High-resolution Optics Module

Abstract: This is a simple tube microscope, using an RMS objective lens and a second lens for tube length correction. It allows the OpenFlexure microscope to take images with higher resolution than those obtained with the webcam lens.

Description: This is the optics module for the OpenFlexure Microscope. It holds a Raspberry Pi camera module and an objective lens. It's compatible with the LS65 or LS75 version of the stage, depending on whether the objective has a 25mm or 35mm parfocal distance.


This is a simple tube microscope, using an RMS objective lens and a second lens for tube length correction. It allows the OpenFlexure microscope to take images with higher resolution than those obtained with the webcam lens.

Requires:

- **Mechanical holder:** [Optics module plastic parts x 1](#)
- **Optics and Sensor:**
- **Tools for removing the lens:** [Camera lens removal tools x 1](#)

Instruction:

Step 1. First, assemble the necessary tools and parts: the optics module plastic part (and camera board cover if you're using a Pi camera), the camera module, the tools for removing the camera lens, the tool for inserting the tube lens, and optionally two M2 screws to secure the camera (if using Pi camera). NB the lens removal tool, board gripper, and optics module are all specific to the camera you're using.

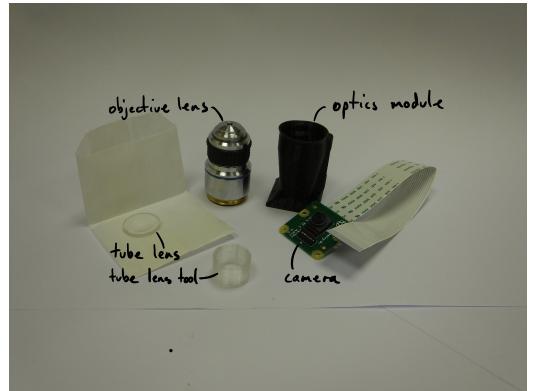


Step 2. WARNING! The camera board is static sensitive. Take the usual anti-static precautions (ideally use an anti-static wristband connected to ground, but at the very least make sure you touch an earthed object, such as a metal pipe, before working on the camera module.




Step 3. Before assembling the parts into the holder, make sure it's free from dust by blowing some air through it, and check there are no strings of plastic in the central hole through the mount.


Step 4. First, insert the tube lens into the holder. Place the lens on the insertion tool, curved side down. Take the optics module, and line it up with the lens. Push down, and the lens should pop into place. Take care to keep the optics module straight as you do this.

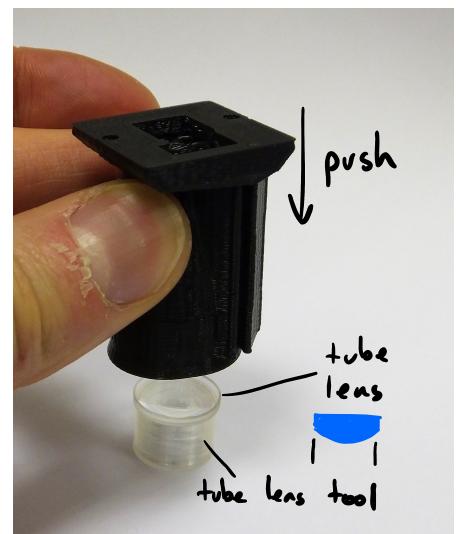



Step 5. If the lens does not sit flat, place the optics module on the insertion tool and push again - it may require some force, and you should hear it click into place.

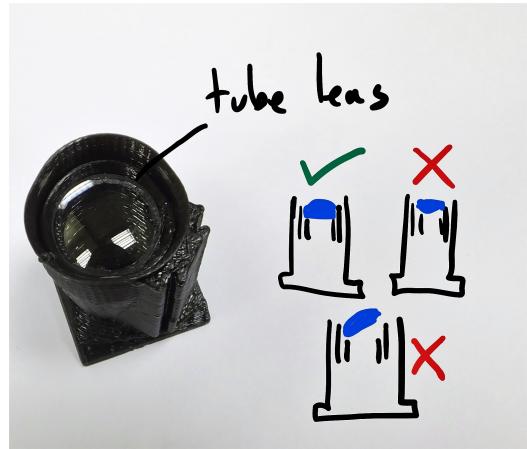

Step 6. Remove the protective film from the camera lens.

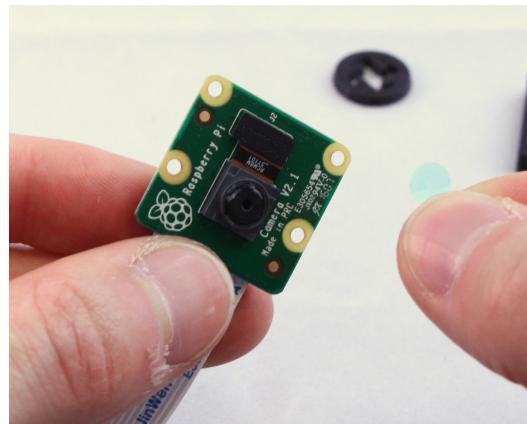


Step 7. There is a small ribbon cable connecting the camera to the PCB that is very easy to break. There is a square plastic jig that fits over the camera and PCB (the "camera board gripper"), which stops the camera twisting and damaging the ribbon cable. Fit this over the camera as shown. Note that the part for v2 of the camera board will sort-of fit v1, but you need to be a little more careful as it's not a perfect fit. For USB webcams, the lens can simply be unscrewed from the back of the PCB with a small screwdriver.

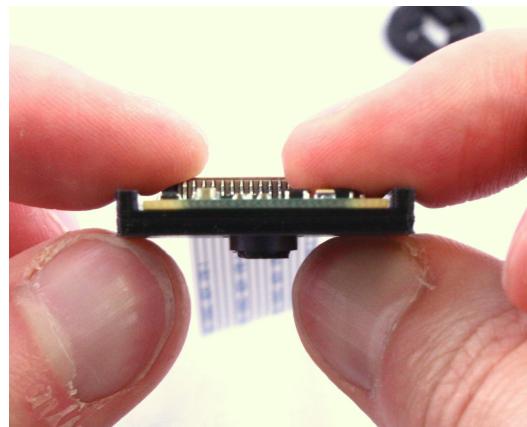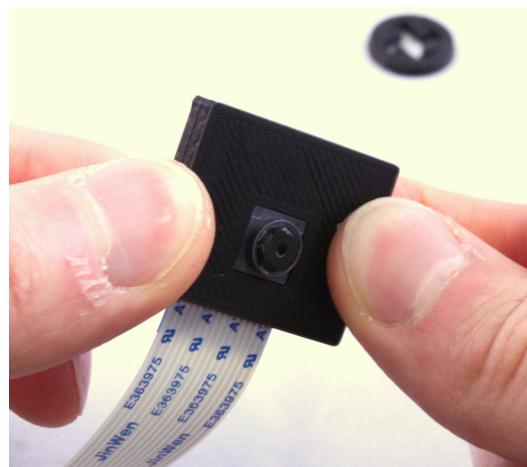


Step 8. Next, unscrew the lens from the camera module. Use the plastic tool to grip the lens module. This is a small circular part with four prongs that fits over the lens of the camera board (version 2 only) as shown. To remove the lens, push the removal tool onto the lens (just the top part, with the little plastic flanges) and turn anticlockwise to remove it. It only works if the prongs are pointing anticlockwise, so make sure it's the right way round. You will may feel a crack as it starts to move - this is normal. It's important to use the board gripper to hold the camera chip in place and prevent damage to the delicate ribbon cable. After you've removed the lens, check that the little black ribbon cable connecting the camera module (square plastic housing) to the PCB is still connected - pop it back in by pushing it with a finger if needed. If you are using a USB camera, the lens is held on by two screws from the back, and can be removed by unscrewing these. **SAVE THE SCREWS**, as you will need them later! Once you've removed the lens, either cover the sensor (without touching it) or fit the camera to the optics module immediately, to avoid dust buildup on the sensor..



Step 9. Press the camera module onto the bottom of the optics mount, so the camera board sits flush with the bottom of the mount - it should fit in with a gentle but firm press from your thumb. It probably won't crack or click into place, but it should be held by the push-fit. Don't worry if it's not gripped very strongly, the cover will hold it in place. If you're using the USB version, the board simply screws onto the optics module using the screws that held the lens on.

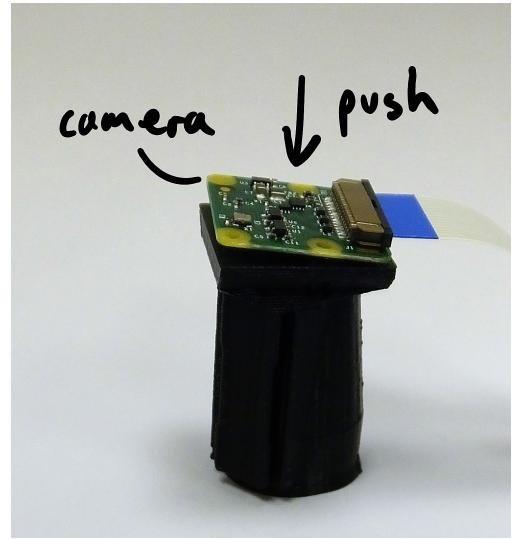

Step 10. Place the cover on top of the camera board, then screw the cover and the camera board onto the optics module, using two M2 screws.

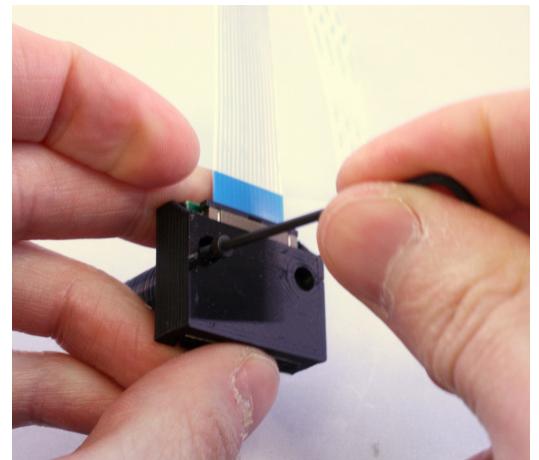

Step 11. Finally, screw the objective into the tapered opening at the top of the module. Take care to keep it aligned with the axis of the module - when it is screwed down, you don't want a gap anywhere between the top of the plastic and the bottom of the objective "shoulder". This may be a tight fit - sometimes clamping the optics module in a vice is helpful.



Actuator assembly

Abstract: An actuator, based on an M3 screw, brass nut, and elastic bands, that can precisely and repeatably move a plastic flexure mechanism.

Description: This assembly is used to actuate each of the three axes of the OpenFlexure Microscope stage. It consists of a plastic gear that fits onto an M3 screw, which then pulls up on a brass nut. Elastic bands pull down on the nut, tensioning the system. This confers two benefits: firstly, it allows us to actuate the flexure mechanism both up and down - doubling the travel - and secondly, it eliminates much of the backlash from the system.

An actuator, based on an M3 screw, brass nut, and elastic bands, that can precisely and repeatably move a plastic flexure mechanism.

License: CERN Open Hardware License

Requires:


- **Lead screw:** 25mm M3 Hexagon-head screw, Qty: 3 x 1
- **Nut for lead screw:** M3 Nuts (preferably brass), Qty: 4 x 1
- **Washer:** M3 Washer, Qty: 3 or 4 x 1
- **Return spring:** Viton O ring, 30mm inner diameter, 2mm cross section, Qty: 3 x 1
- **Actuating gear/thumbwheel:** Gear/thumbscREW, Qty: 3 x 1
- **Microscope Foot/anchor for return spring:** Microscope feet, Qty: 3 x 1

Instruction:


Step 1. There are three of these actuators on the microscope: these steps should be followed for each one in turn. Note the quantities in the parts list are for one actuator, and they don't include the microscope body. The image for this step shows everything you will need to assemble the three actuators (NB the illumination arm is not yet needed, and the band insertion tool is now longer and double-ended).




Step 2. We need to place a nut inside each actuator. Insert it from the side, through the slot in the actuator column, and push it into place with the tool. Take care to insert the nut corner-first (i.e. with flat sides parallel to the tool), otherwise it will jam.

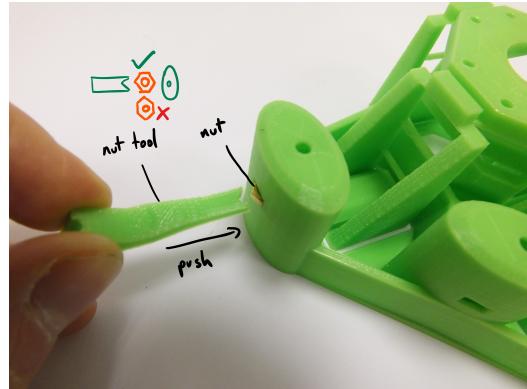

Step 3. Next, put one of the hex-head screws into the gear, and put a washer on the screw thread underneath the gear. Put a very small blob of grease (if available) on the screw thread. Insert the screw into the hole at the top of the column, and screw it into the nut that is now inside the actuator. Use a washer between the gear and the microscope body.


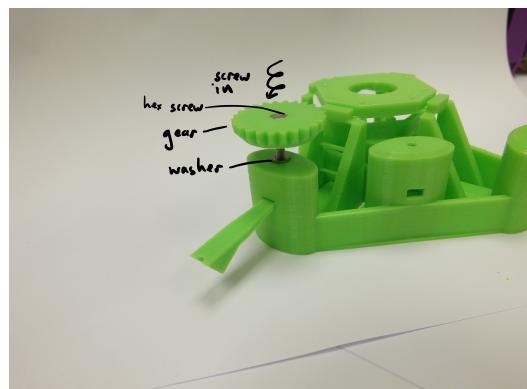


Step 4. Turn the microscope over and, using the elastic band tool or a craft knife, snap the two thin plastic supports as shown if they are present.


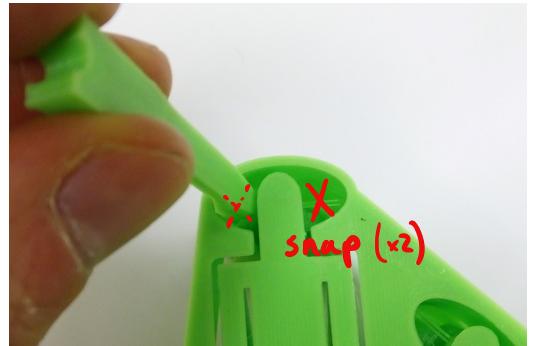


Step 5. Fit the foot to the bottom of the microscope. NB the two tilted feet must go on the outside of the microscope, and the gap in the foot points towards the optics module.

Step 6. Hook one of the Viton (or rubber) bands over the elastic band tool, and insert it through the foot. You should end up with one loop of elastic band passing over each fork-shaped end of the insertion tool. The flat side of the tool should be on the outside, so that the channels in either end face towards the centre of the foot.

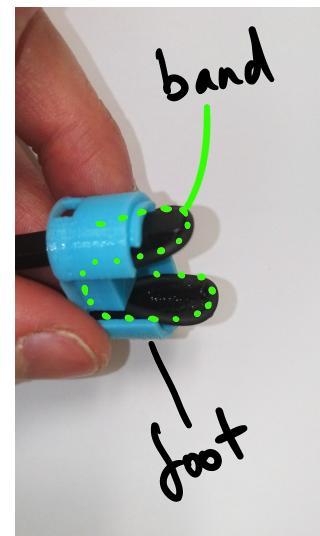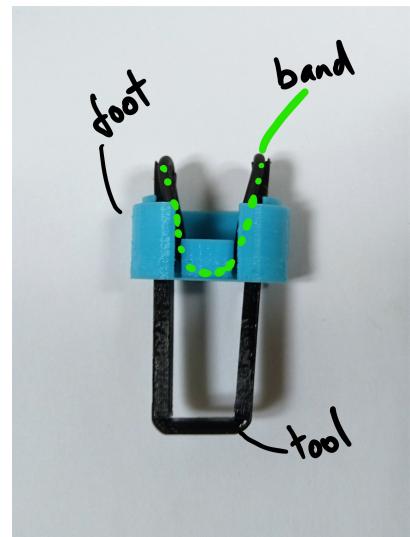

Step 7. Place the foot on the body, taking care that it's in the right place (the tilted ones are for the outer two actuators, the flat one is for the centre) and the right way round (the open part of the foot points towards the optics module in the centre of the microscope). Push the insertion tool in until the band clicks into place. You may find it helps to push on the gear or use the nut insertion tool to hold the actuator in place while you do this. What's happening here is that the elastic band is pushed over little hooks inside the actuator column, near the top, and the click is the band popping into place. The band pulls down on the actuator, so after it's in place the gear should turn with a little more resistance - this is important for the microscope to work properly.

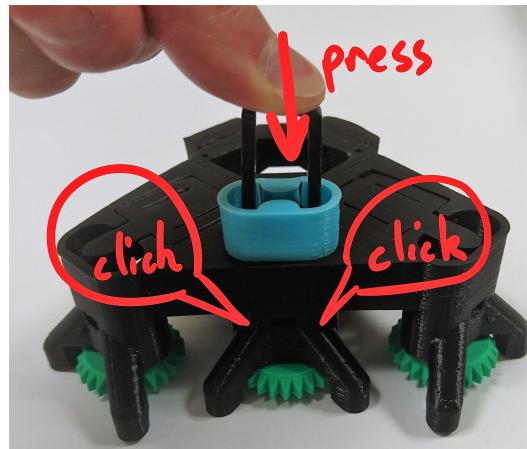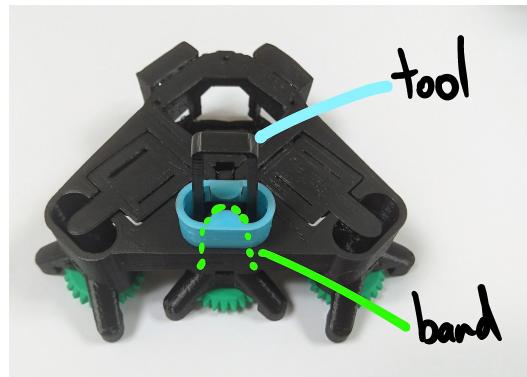

Step 8. You should now have the band pulling both sides of the actuator down as shown. If you used a long elastic band, you can break the thin layer at the bottom of the foot on either side, to pull the band out and wrap it around the foot to tighten it. This is not necessary for the Viton bands, or if you're using doubled-over #32 bands.



Step 9. Congratulations, you've assembled the actuator! There are three to do, so you may now need to go back and make the next one.

Illumination Arm (fixed)

Abstract: A holder for an LED to take images in transmitted light.

Description: The illumination arm holds an LED above the sample, for simple transmission images. The adjustable version with a condenser lens will get higher resolution and brightness.



A holder for an LED to take images in transmitted light.




Requires:



- **Plastic mounting parts:** [Illumination arm and rear foot x 1](#)
- **Wiring:** [Jumper wires with female header pin connectors x 1](#)
- **LED:** [White LED, 5mm diameter, 15 degree beam angle x 1](#)
- **Attachment screw:** [8mm M3 screws, Qty: 3 x 1](#)
- **Washer for illumination arm:** [M3 Washer, Qty: 3 or 4 x 1](#)
- **Nut for illumination:** [M3 Nuts \(preferably brass\), Qty: 4 x 1](#)

Instruction:

Step 1. Assemble the illumination: Bend the LED so it's pointing at right angles to the cable, then thread the wire for the LED through the illumination arm and push the LED into its hole. It's usually easiest to start at the LED end and feed the connector into the mount (it might require a little force). If you have a pre-made LED wire with a connector on the end, threading it through is all that's needed - if you just have the components you will need to connect them together in series first, and solder or crimp them. Some of the pre-made LED wires are thicker than expected, so they take some force to pull into the arm - but they will go in with some perseverance! You can stick a piece of white tape or sticky label over the LED (i.e. between the LED and the stage) to act as a diffuser, this helps to blur out out-of-focus artefacts, such as dirt on the wrong side of your microscope slide.

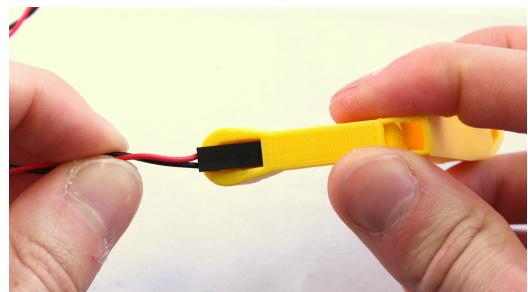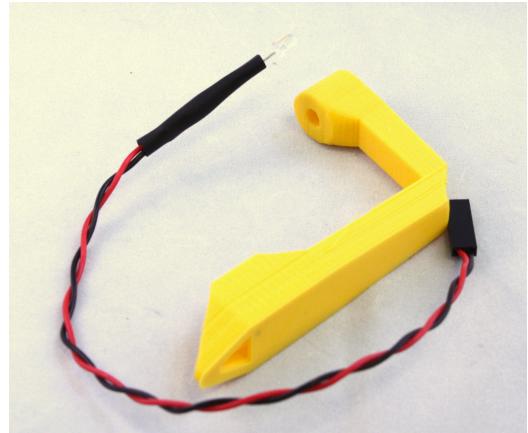

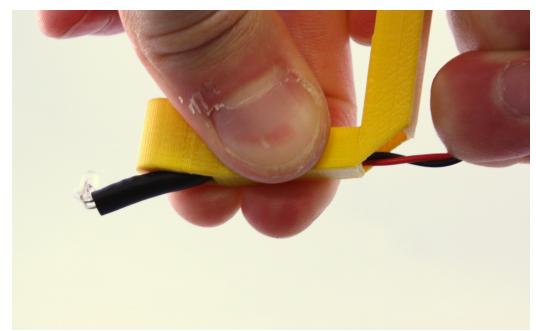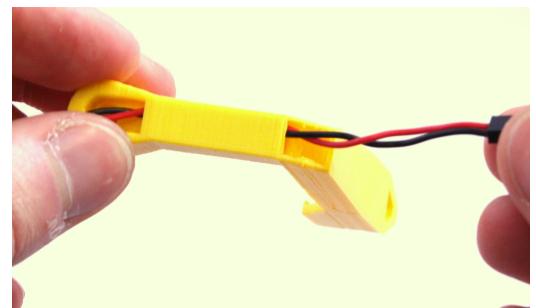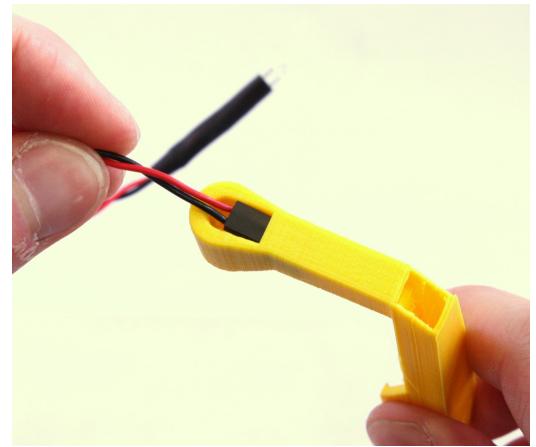

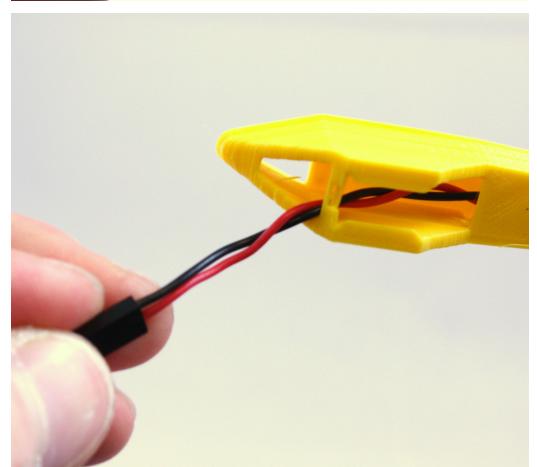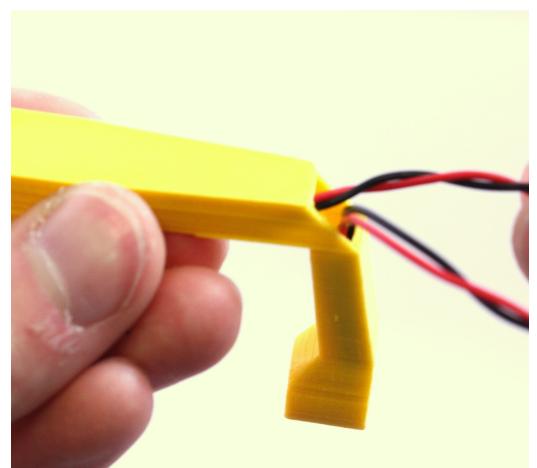

Step 2. Slot the illumination arm on to the dovetail mount for it at the back of the microscope. Wiggling it slightly makes it easier to push on. It should be pushed on from the bottom, and on version 5.15.1 and later, there should be a stop to prevent it sliding up and off.



Adjustable Illumination Arm

Abstract: A holder for an LED to take images in transmitted light.

Description: The illumination arm holds an LED above the sample. There is a screw holding the two parts together, which goes through a slot to allow simple, slip plate style alignment of the LED. You can use some white tape stuck over the LED as a diffuser to improve the microscope's resolution - though you may want to consider using the condenser lens version if you have an RMS microscope objective.

A holder for an LED to take images in transmitted light.

Requires:


- **Plastic mounting parts:** Illumination arm and rear foot x 1
- **Wiring:** Jumper wires with female header pin connectors x 1
- **LED:** White LED, 5mm diameter, 15 degree beam angle x 1
- **Attachment screw:** 8mm M3 screws, Qty: 3 x 1
- **Washer for illumination arm:** M3 Washer, Qty: 3 or 4 x 1
- **Nut for illumination:** M3 Nuts (preferably brass), Qty: 4 x 1

Instruction:

Step 1. To build the condenser, you will need the two parts of the mounting arm, an M3 screw and washer, the cable, and the LED (ignore the condenser lens shown in the picture).




Step 2. Assemble the cable for the LED. If using pre-crimped cables, as supplied with our kit, fit the 2-way header to either end. Don't worry about polarity, it's easy and harmless to simply rotate the connector 180 degrees if you need to.

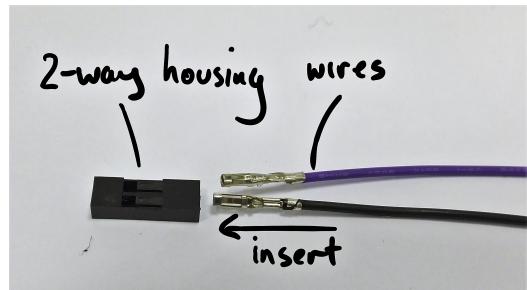




Step 3. Thread the cable through the larger, vertical part of the mounting arm as shown. The smaller part is different from the one shown in the picture, but the cable should still fit through.




Step 4. Screw the horizontal part of the mounting arm onto the vertical part. You should be able to put a nut into the top of the vertical part.

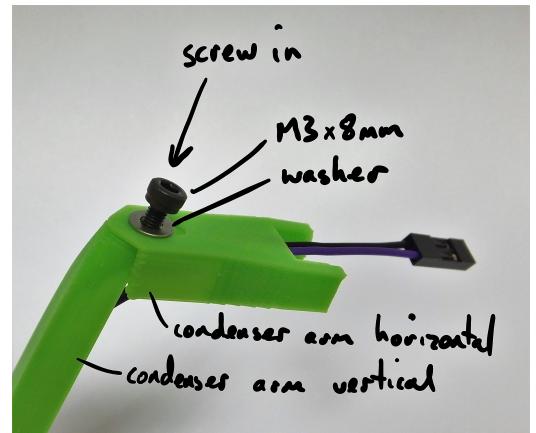

Step 5. Connect the LED to the cable by pushing the legs into the 2-way housing. Bend the lens part of the LED at right angles to the legs, lift it up, and insert it into the hole.

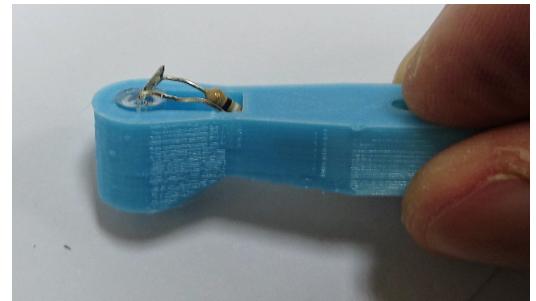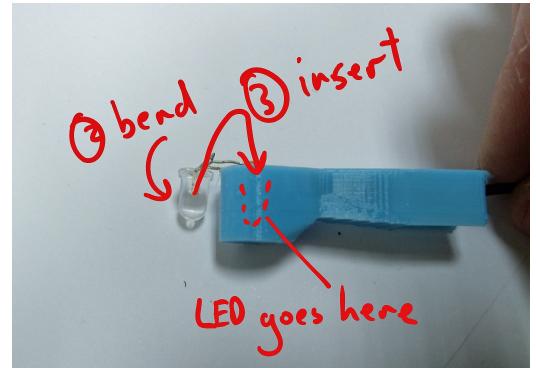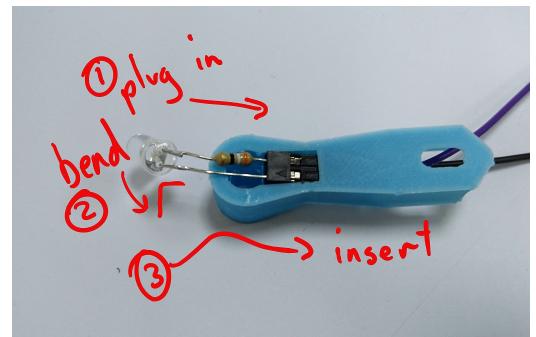

Step 6. Fit the illumination arm to the microscope by sliding it onto the dovetail from underneath. The thin linker between the two sides of the clip may snap - this is not a problem, it's there only for support during printing..

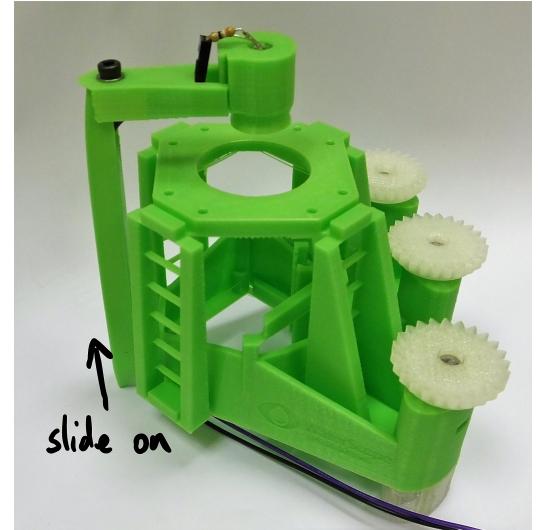

Illumination Arm with Condenser

Abstract: A holder for an LED and lens to take images in transmitted light.

Description: The illumination arm holds an LED above the sample, with a condenser lens. It will get higher resolution and brightness. There is a screw holding the two parts together, which goes through a slot to allow simple, slip plate style adjustment of the condenser alignment. There is also a simpler fixed arm holding a single LED.

A holder for an LED and lens to take images in transmitted light.

Requires:


- **Plastic mounting parts:** Illumination arm and rear foot x 1
- **Wiring:** Jumper wires with female header pin connectors x 1
- **LED:** White LED, 5mm diameter, 15 degree beam angle x 1
- **Attachment screw:** 8mm M3 screws, Qty: 3 x 1
- **Washer for illumination arm:** M3 Washer, Qty: 3 or 4 x 1
- **Nut for illumination:** M3 Nuts (preferably brass), Qty: 4 x 1
- **Condenser lens:** Condenser Lens (optional) x 1

Instruction:


Step 1. To build the condenser, you will need the two parts of the mounting arm, an M3 screw and washer, the cable, LED, and lens, and the condenser housing. You will also need the printed lens insertion tool to put the lens into the holder.





Step 2. Place the condenser lens on the insertion tool (flat side down), and push the condenser housing down onto it. Take care to keep it vertical. You may need to push fairly hard. As with the tube lens in the optics module, check it's flat and push again if necessary.

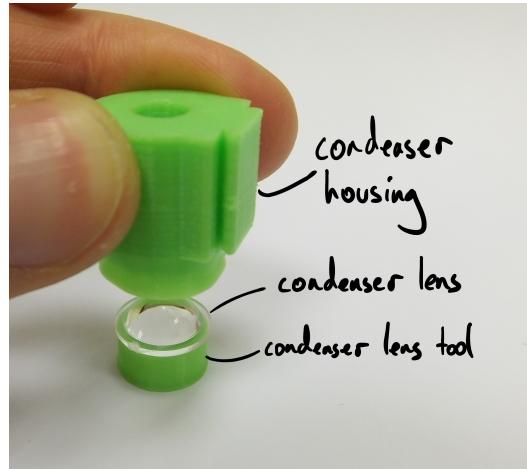

Step 3. Insert the LED into the condenser housing. It should push-fit and stay in place. It can be glued if it doesn't stay in.

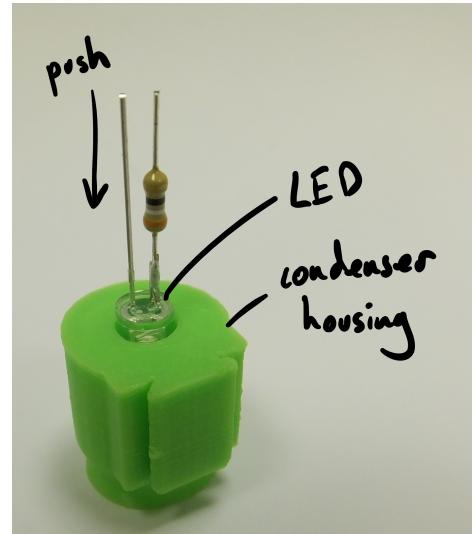


Step 4. Assemble the cable for the LED. If using pre-crimped cables, as supplied with our kit, fit the 2-way header to either end. Don't worry about polarity, it's easy and harmless to simply rotate the connector 180 degrees if you need to.

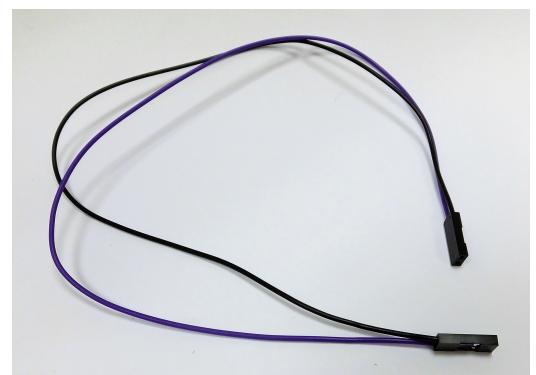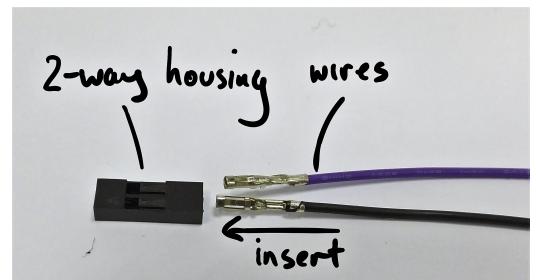


Step 5. Thread the cable through the mounting arm as shown. It's much easier to do this before assembling the arm.



Step 6. Screw the horizontal part of the mounting arm onto the vertical part. The hole should be about the right size for an M3 screw to self-tap.

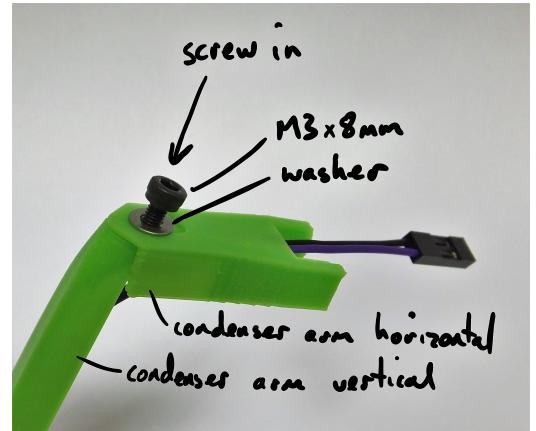


Step 7. Connect the LED to the cable by pushing the legs into the 2-way housing.

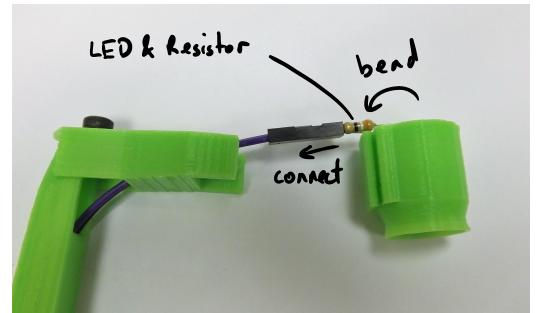

Step 8. Slide the condenser housing onto the dovetail clip at the end of the arm. Wiggling it from side to side might help if it's stiff.



Step 9. Fit the illumination arm to the microscope by sliding it onto the dovetail from underneath. The thin linker between the two sides of the clip may snap - this is not a problem, it's there only for support during printing..

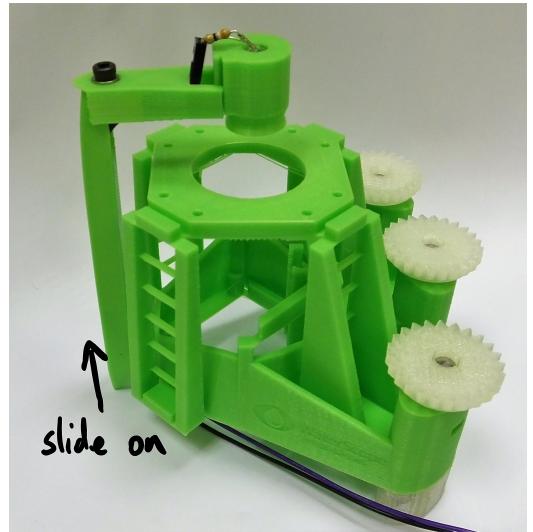

Part: M3 Nuts (preferably brass), Qty: 4

Description: M3 Nuts (preferably brass), Qty: 4

Supplier: Anglian Fasteners Limited

Supplier catalog #: 11355

URL: www.anglianfasteners.co.uk

Instruction:

Step 1. These nuts form part of the actuator. Brass ones will have less friction and result in smoother motion. An extra one (doesn't have to be brass, but it's fine if it is) is required by the adjustable illumination arm.

Part: 25mm M3 Hexagon-head screw, Qty:3

Description: 25mm M3 Hexagon-head screw, Qty:3

Supplier: Anglian Fasteners Limited

Supplier catalog #: 6663

URL: www.anglianfasteners.co.uk

Instruction:

Step 1. These screws are used as the lead screws for the stage; stainless steel is probably best. If hexagon-head screws are not available, a good substitute is to screw a nut tightly against the head of the screw, so the screw can still be turned by the thumbwheel or gear.

Part: M3 Washer, Qty: 3 or 4

Description: M3 Washer, Qty: 3 or 4

Supplier: Anglian Fasteners Limited

Supplier catalog #: 7118

URL: www.anglianfasteners.co.uk

Instruction:

Step 1. The washers fit between the gears and the microscope. An extra one is required by the adjustable illumination arm.

Part: 8mm M3 screws, Qty: 3

Description: 8mm M3 screws, Qty: 3

Supplier: Anglian Fasteners Limited

Supplier catalog #: 14315

URL: www.anglianfasteners.co.uk

Instruction:

Step 1. These screws attach the sample clips. A further screw is used in the illumination arm, if you have an adjustable one.

Part: White LED, 5mm diameter, 15 degree beam angle

Description: White LED, 5mm diameter, 15 degree beam angle

Supplier: Farnell Element14

Supplier catalog #: 1716696

Manufacturer catalog #: OVL-3321

Instruction:

Step 1. The white LED is for illumination. Exact specifications are unimportant, but this is the one I use. 5mm diameter is what I've designed for, but there's no reason you couldn't tweak the illumination arm to use a 3mm LED.

Part: 60 Ohm resistor

Description: 60 Ohm resistor**Supplier:** Farnell Element14**Supplier catalog #:** 1565328**Manufacturer catalog #:** MFR3-68RFC

Instruction:

Step 1. This resistor allows the LED to be run from a 5V GPIO pin on the Raspberry Pi. If you use a different LED, you may need a different resistor. Note the maximum current you can draw from the Pi is not huge.

Part: Jumper wires with female header pin connectors

Description: Jumper wires with female header pin connectors

Instruction:

Step 1. You need wires to connect the LED to the 5V power supply from the Raspberry Pi. If you've bought a kit, it will be supplied with a ready-made lead that has the LED at one end and the header connector at the other end. My preferred DIY option is two cores from a ribbon cable, soldered to the LED and resistor at one end and with a crimped-on connector at the other. If you want a no-soldering solution, four female-female jumpers allow you to connect everything (using a cut-off leg of the resistor to join two of them together).

Part: Viton O ring, 30mm inner diameter, 2mm cross section, Qty: 3

Description: Viton O ring, 30mm inner diameter, 2mm cross section, Qty: 3

Instruction:

Step 1. The O-rings function as springs, to provide a return force on the nut that is pulled upwards by the screw. I use Viton O-rings, 30mm ID, as they provide about the right force. You can substitute for number 32 elastic bands, which are about 3mm wide and about 150mm in circumference (i.e. 75mm long when measured as a loop). If you use these, you should fold them in half, so you have a band that's effectively 75mm in circumference. Buying a band that size is a bit difficult, but if you can get one, great!

Part: Raspberry Pi

Description: Raspberry Pi

Instruction:

Step 1. Any model will do! If you buy a kit, it may or may not come with a Pi. NB if you buy a Pi Zero, only the Zero+ has a camera connector, and even that needs a special lead. Of course, if you're using a camera that's not the Raspberry Pi camera module, you can use whatever USB-compatible device you like.

Part: Camera Module

Description: Camera Module



Supplier: Raspberry Pi

Manufacturer catalog #: Camera Module v2

Instruction:

Step 1. NB I've pictured v1.3 of the camera board here; the microscope is designed to work with v2 of the board, although older versions of the optics module exist that work with v1. You can also use a webcam, and there are versions of the optics module available to use the Logitech C270 and a USB webcam supplied by WaterScope.

Part: Microscope Body

Description: Microscope Body

File: [main_body_SS.stl](#)

File: [main_body_SS-M.stl](#)


File: [main_body_LS65.stl](#)

File: [main_body_LS65-M.stl](#)

Material usage: 0.1 KG

Instruction:


Step 1. This should be printed without support material. On smaller/less well calibrated machines, I print this part on its own and then print the rest of the parts in a second print. Which version of the body you need depends on two things: firstly, whether you want to use a large optics module (versions with LS65 in the name) or a small one (versions with SS), and secondly whether you want attachment lugs for motors (files with -M in the name). All the files start with 'body_'. You need the SS version if you're using a Raspberry Pi camera or a Logitech C270 camera, together with the lens that came on the camera. If you're using a microscope objective, or the USB camera with M12 lens, you need the LS65 version. The microscope body takes around 8 hours on a RepRapPro Ormerod (and many other low-end printers) or about 5 hours on Ultimaker, MakerBot, and the like.

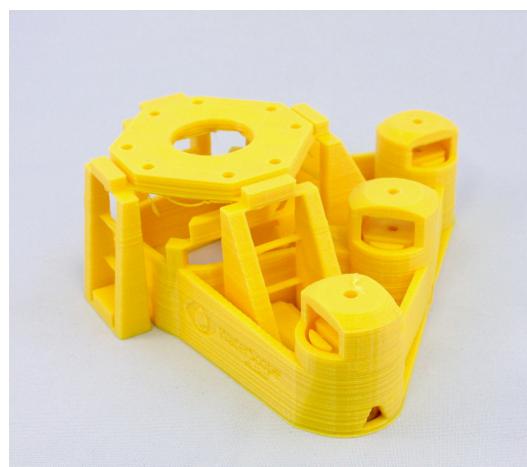

Step 2. After printing, you should run a 3mm drill bit through the 3mm holes in each actuator, to ensure the screws can rotate freely. If you don't have a drill, an M3 screw should do...



Step 3. If the bottom layer has oozed out too much, or if you used a brim (not recommended) you might need to clean up the bottom so it looks like this - it's important that the moving parts aren't stuck to the body with a thin layer of plastic (this most often happens around the struts connecting the objective clip to the rest of the microscope).

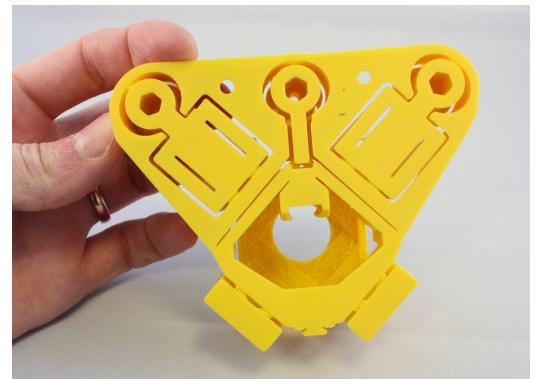

Step 4. You may need to use needle-nose pliers to pull strings of plastic from the underside of the microscope stage or the underside of the caps of the actuator columns. If your printer is correctly calibrated there shouldn't be much, and I often get away without any - but some printers (particularly if using ABS) are prone to a bit of "spaghetti" under the stage.

Part: Microscope feet, Qty: 3

Description: Microscope feet, Qty: 3

File: [feet.stl](#)

Material usage: 0.01 KG

Instruction:

Step 1. You need two tilted feet and one untilted foot, as provided in feet.stl.

Part: Optics module plastic parts

Description: Optics module plastic parts

File: [optics_picam2_pilens.stl](#)

File: [optics_picam2_rms.stl](#)

File: [picam_cover.stl](#)

File: [optics_usbcam_M12.stl](#)

File: [optics_usbcam_rms.stl](#)

File: [optics_c270_owlens.stl](#)

[**File: optics_c270_rms.stl**](#)**Material usage:** 0.01 KG

Instruction:

Step 1. There are several versions of the optics module, depending on your camera (Raspberry Pi Camera v2, Logitech C270, or WaterScope USB camera) and on whether you will use the lens from the camera (pilens, M12, ownlens) or an RMS objective and 40mm tube lens. Make sure you pick the right STL file for your camera module! There is a cover that fits over the Raspberry Pi camera module, and holds it firmly onto the optics module.

Step 2. The optics module needs to print with some fine detail, so the dovetail meshes nicely with the stage. A good way to ensure this is to print it at the same time as other parts - either print more than one optics module at a time, or print it at the same time as the microscope body. This slows down the time for each layer, and means the plastic can cool more completely before the layer on top is deposited, resulting in a higher-quality part. The optics module is best printed in black to cut down on stray light inside the tube - though it will still work in other colours.

Part: Camera lens removal tools

Description: Camera lens removal tools[**File: picam2_lens_remover.stl**](#)[**File: picam2_board_gripper.stl**](#)[**File: picam1_board_gripper.stl**](#)**Material usage:** 0.01 KG

Instruction:

Step 1. These tools make it much easier to remove the lens from a Raspberry Pi camera module. See the documentation for the Raspberry Pi Camera Mount/Extension Tube for instructions. These tools are for version 2 of the camera board. The board gripper for v1 is included for convenience - you'll need a small pair of pliers (I like mini combination pliers) to unscrew the lens if you have v1 of the board. v2 of the Raspberry Pi camera module now comes with a lens removal tool, which is better than the printed one - but the board gripper is still useful.

Part: Gear/thumbscrew, Qty: 3

Description: Gear/thumbscrew, Qty: 3[**File: large_gear.stl**](#)**Material usage:** 0.01 KG

Instruction:

Step 1. These gears are used to control the leadscrews that move the stage. You can replace them with more comfortable thumbscrews if desired, or use stepper motors to drive them automatically. You need three gears. NB the gears are designed to print with the circular face down, this avoids marring of the teeth and helps it stick to the bed better.

Part: Optics (optional)

Description: Optics (optional)

Instruction:

Step 1. If you are not simply using the lens that came with your camera, you will need an RMS microscope objective. The LS65 version of the microscope body is compatible with 25mm parfocal length objectives; if you have a longer 35mm parfocal length objective you'll need to either print a taller LS75 body, or fit a 10mm spacer to raise the sample up higher. If you are using an RMS objective, you'll also need to use a tube lens; 16mm diameter, 40mm focal length (e.g. [Comar 40 PQ 16] (http://www.comaroptics.com/components/lenses/simple-convex-lenses/quality-planoconvex-lenses#row-40_pq_16))

Part: Condenser Lens (optional)

Description: Condenser Lens (optional)

Instruction:

Step 1. If you want to use a condenser lens, so far I've based the design around a 12mm diameter lens, with a focal length of about 12mm. WaterScope can supply these, contact them for more details.

Part: Illumination arm and rear foot

Description: Illumination arm and rear foot

File: [illumination_and_rear_foot_LS65.stl](#)

File: [illumination_and_rear_foot_LS65_condenser.stl](#)

File: [illumination_and_rear_foot_LS75_condenser.stl](#)

File: [illumination_and_rear_foot_LS75_condenser_tall.stl](#)

File: [illumination_and_rear_foot_SS40.stl](#)

Material usage: 0.01 KG

Instruction:

Step 1.

The illumination arm fits on to the microscope using a dovetail, and includes the rear foot. The condenser lens version will also work without a lens, though the image won't be as bright.

The illumination arm is available in several versions:

`illumination_and_rear_foot_[version].stl`. You can choose either adjustable arm with a single LED or one that uses a condenser lens (`_condenser`). It's important to match the type of microscope you're using (`LS65`, `LS75`, or `SS40`). If you are using a riser to make space for a larger objective (if you've got a 35mm parfocal length objective) you'll need the `LS75` version. If you are using tall feet, you will need the `_tall` version.

Part: Sample Clips

Description: Sample Clips

File: [sample_clip.stl](#)

Material usage: 0.002 KG

Instruction:

Step 1. These are optional, but they are useful for holding slides on the stage. Usually, I use two of these for a normal microscope slide, screwed onto the top of the stage.

Part: Microscope leg test object

Description: Microscope leg test object

File: just_leg_test.stl

Material usage: 0.01 KG

Instruction:

Step 1. This test piece checks your printer can bridge between the tops of the legs without failing. If it passes this test, you're probably fine to print the rest of the microscope.

Authors

| Name | E-mail | Affiliation | ORCID |
|-----------------------|---------------------------|-------------------------|---------------------|
| Dr. Richard W. Bowman | richard.bowman@cantab.net | University of Cambridge | 0000-0002-1531-8199 |

Total bill of materials for this project

| Part | Quantity | Supplier | Supplier part number | URL |
|---|----------|---------------------------|----------------------|--|
| <u>Optics module plastic parts</u> | 2 | | | |
| <u>Camera lens removal tools</u> | 2 | | | |
| <u>Illumination arm and rear foot</u> | 3 | | | |
| <u>Jumper wires with female header pin connectors</u> | 3 | | | |
| <u>White LED, 5mm diameter, 15 degree beam angle</u> | 3 | Farnell Element14 | 1716696 | |
| <u>8mm M3 screws, Qty: 3</u> | 4 | Anglian Fasteners Limited | 14315 | www.anglianfasteners.co.uk |
| <u>M3 Washer, Qty: 3 or 4</u> | 6 | Anglian Fasteners Limited | 7118 | www.anglianfasteners.co.uk |
| <u>M3 Nuts (preferably brass), Qty: 4</u> | 6 | Anglian Fasteners Limited | 11355 | www.anglianfasteners.co.uk |
| <u>Condenser Lens (optional)</u> | 1 | | | |
| <u>25mm M3 Hexagon-head screw, Qty:3</u> | 3 | Anglian Fasteners Limited | 6663 | www.anglianfasteners.co.uk |
| <u>Viton O ring, 30mm inner diameter, 2mm cross section, Qty: 3</u> | 3 | | | |
| <u>Gear/thumbscrew, Qty: 3</u> | 3 | | | |
| <u>Microscope feet, Qty: 3</u> | 3 | | | |
| <u>Microscope Body</u> | 1 | | | |
| <u>Sample Clips</u> | 1 | | | |
| <u>Raspberry Pi</u> | 1 | | | |
| <u>Microscope leg test object</u> | 1 | | | |

